

SeqTrimNext

Statistics of pre-processing

Plataforma Andaluza de Bioinformática
Universidad de Málaga

April 24, 2020

1 Output Files

SeqTrimNext provides several files, the most interesting ones are in the following directories:

- `output_files`
 - `output.less`, containing an extensive information about the trimming of each sequence. It can be visualised on terminal using the command `less -R`.
 - `used_params.txt`, containing the complete set of parameters used for execution of SeqTrimNext with your data
 - `rejected.txt`, containing a list of rejected sequences together with the reason for their removal.
 - `initial_stats.json`, containing statistics for raw sequences.
 - `stats.json`, containing the statistics of the cleaning process.
 - There is a collection of `folders` that gather sequences with the same MID; each folder contains a `sequences` file (in FASTQ format) with useful reads. There may also exist a file with reads containing low complexity regions. If you want to reconstruct a SFF with the useful segment of each pre-processed read, use `sff_info` file in combination with the original SFF file for the `sfffile` tool.
- `graphs`
 - `size_stats.png`, a graph with the distribution of read lengths in raw data (see Fig. ??).
 - `qualities.png`, a graph to inspect read qualities in raw data (see Fig. ??).
 - `PluginExtractInserts_insert_size.png`, a graph with the distribution of read lengths after SeqTrimNext pre-processing (see Fig. 1).
 - There are other graphs (mostly bar plots) that illustrate the quality of pre-processed reads. All are in PNG format.
- `latex`
 - It is provided as a compressed file `latex.zip` containing all “.tex” files required to compile this document. Graphs are taken from the `graph` folder

2 Relevant parameters

In this section, the relevant parameters used in your experiment are shown. Full information about the parameters can be obtained from file `used_params.txt`

2.1 General

Plugins applied to every sequence, separated by commas. Order is important

1. PluginIndeterminations
2. PluginFindPolyAt
3. PluginAbAdapters
4. PluginUserContaminants
5. PluginContaminants
6. PluginVectors
7. PluginLowQuality
8. PluginLowComplexity
9. PluginExtractInserts

Remove duplicated (clonal) sequences (using CD-HIT 454)

```
remove_clonality: false
```

Minimum insert size for every trimmed sequence

```
min_insert_size_trimmed: 30
```

Minimum insert size for each end of paired-end reads; true paired-ends have both single-ends longer than this value

```
min_insert_size_paired: 40
```

Seqtrim version

```
seqtrim_version: 2.0.67
```

```
min_sequence_size_raw:
```

2.2 Quality

Minimum quality value for every nucleotide

```
min_quality: 20
```

```
window_width:
```

2.3 Contaminants

Blast E-value used as cut-off when searching for contaminations

```
blast_evalue_contaminants: 1.0e-10
```

Minimum required identity (%) for a reliable contamination

```
blast_percent_contaminants: 85
```

Minimum hit size (nt) for considering a true contamination

```
min_contam_seq_presence: 40
```

Genus of input data: contaminations belonging to this genus will be ignored

```
genus:
```

Is a contamination considered a source of sequence rejection? (setting to false will only trim contaminated sequences instead of rejecting the complete read)

```
contaminants_reject: true
Path for contaminants database
    contaminants.fasta
    cont_ribosome.fasta
```

3 Pre-processing statistics

Next figure is equivalent to Figure ?? but using output reads (useful sequences). The mode is expected to decrease but the shape of the plot should be similar.

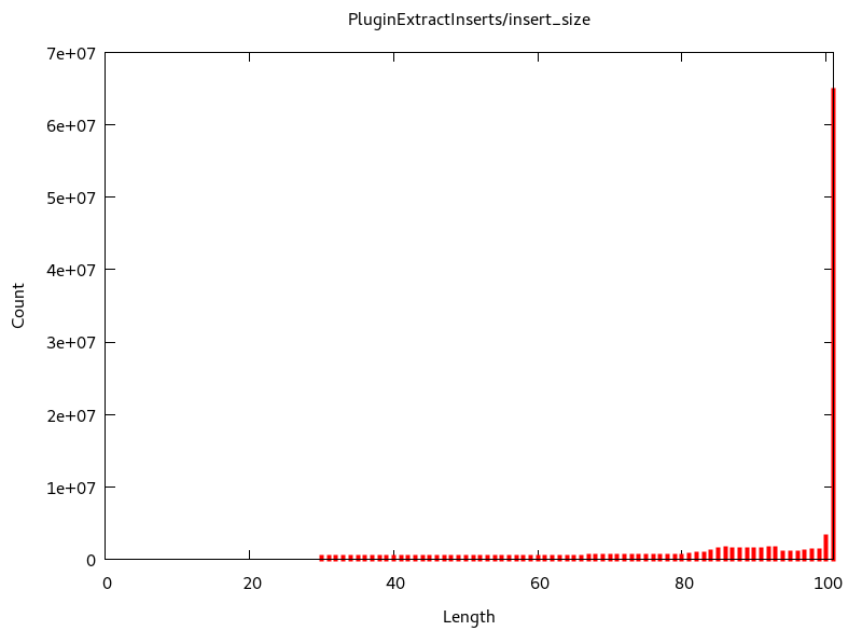


Figure 1: Size distribution of the output sequences. Short sequences ($< \text{min_insert_size_trimmed}$) were removed. [PluginExtractInserts_insert_size.png]

Summary statistics of the SeqTrimNext analysis. Be careful and read all warnings that are indicating concerns about your data. In the files `initial_stats.json` and `stats.json` can be found a full statistics of your data and SeqTrimNext pre-processing

Input reads:	total	121440546
	Smallest read (bp)	101
	Largest read (bp)	101
	Mode (bp)	0
	Mean (bp)	0.0
Output results:	total	4896067
	Rejected	14103440
	Low complexity reads	82327
	Mode (bp)	91
	Mean (bp)	91.6
	Output paired reads	102358712
	Total output reads	107254779
Linkers:		

Table 1: List of the most frequent Vectors found among your reads

Vectors	sequences
Cloning vector pAAV-MCS, complete sequence.	593258
Cloning vector pKOHPRT complete sequence.	583649
Enterobacteria phage lambda	176591
Cloning vector pVLH/hsp	54165
Cloning vector pWormgate2, complete sequence.	50743

Table 2: List of the most frequent Adapters found among your reads

Adapters	sequences
ABI_Solid3_Adapter_A	344351
ABI_Solid3_GAPDH_Reverse_Primer	289327
TruSeq_Universal_Adapter	281525
ABI_Solid3_GAPDH_Forward_Primer	252216
Illumina_Single_End_Adapter_1	249430

Table 3: List of the most frequent Contaminants found among your reads

Contaminants	sequences
rRNA_small_subunit_Metazoa_Myotis	162309
Saccharomyces cerevisiae S288c chromosome XII, complete sequence	125918
rRNA_long_subunit_Metazoa_Dasytus	122488
Schizosaccharomyces pombe 972h- chromosome III, complete sequence	94507
Podospora anserina S mat+ unordered scaffolds, whole genome shotgun sequence	93200

Table 4: Summary of nucleotides removed in every plugin.

Plugin	Nucleotides	Percent	Warnings
Low Quality	1711085947	Inf %	OK
Low Complexity	89407405	Inf %	OK
Poly T	35887826	Inf %	OK
Poly A	41425496	Inf %	OK
Contaminants	89873938	Inf %	OK
Adapters	111681141	Inf %	OK
Vectors	68877216	Inf %	OK
Indeterminations	69547	Inf %	OK
Inserts	9834457753	Inf %	iW1

iW1 Warning!, only Inf % of nucleotides are useful

4 Rejected reads

Input sequences	121440546
Output sequences	4896067
Rejected sequences	14103440
Output paired sequences	102358712
Total output sequences	107254779
Low complexity sequences	82327

Table 5: Summary of reads removed in every plugin.

Case	Number of sequences	Percent	Warnings
	14103440	11.613 %	OK
Short inserts	9051689	7.454 %	OK
Empty Inserts	4063982	3.346 %	rdW3
Contaminants	960390	0.791 %	rdW4
No Valid Inserts	19243	0.016 %	OK
Low Complexity	7585	0.006 %	OK
Indeterminations	436	0.000 %	OK
Unexpected Vector	115	0.000 %	OK
Total rejected	14103440	11.613 %	OK

rdW3 Warning!, a 3.346 % of your sequences are empty (without an insert)

rdW4 Warning!, a 0.791 % of your sequences are from a contaminant organism or from organelles

References

- [1] Falgueras et al. SeqTrim: a high-throughput pipeline for preprocessing any type of sequence reads. *BMC Bioinformatics* 11:38 (2010)
- [2] Weizhong Li & Adam Godzik. Cd-hit: a fast program for clustering and comparing large sets of protein or nucleotide sequences. *Bioinformatics* (2006) 22:1658-9

Thanks you for use SeqTrimNext! Send us any comment to scbi support