

SeqTrimNext

Statistics of pre-processing

Plataforma Andaluza de Bioinformática
Universidad de Málaga

April 23, 2020

1 Output Files

SeqTrimNext provides several files, the most interesting ones are in the following directories:

- `output_files`
 - `output.less`, containing an extensive information about the trimming of each sequence. It can be visualised on terminal using the command `less -R`.
 - `used_params.txt`, containing the complete set of parameters used for execution of SeqTrimNext with your data
 - `rejected.txt`, containing a list of rejected sequences together with the reason for their removal.
 - `initial_stats.json`, containing statistics for raw sequences.
 - `stats.json`, containing the statistics of the cleaning process.
 - There is a collection of `folders` that gather sequences with the same MID; each folder contains a `sequences` file (in FASTQ format) with useful reads. There may also exist a file with reads containing low complexity regions. If you want to reconstruct a SFF with the useful segment of each pre-processed read, use `sff_info` file in combination with the original SFF file for the `sfffile` tool.
- `graphs`
 - `size_stats.png`, a graph with the distribution of read lengths in raw data (see Fig. ??).
 - `qualities.png`, a graph to inspect read qualities in raw data (see Fig. ??).
 - `PluginExtractInserts_insert_size.png`, a graph with the distribution of read lengths after SeqTrimNext pre-processing (see Fig. 1).
 - There are other graphs (mostly bar plots) that illustrate the quality of pre-processed reads. All are in PNG format.
- `latex`
 - It is provided as a compressed file `latex.zip` containing all “.tex” files required to compile this document. Graphs are taken from the `graph` folder

2 Relevant parameters

In this section, the relevant parameters used in your experiment are shown. Full information about the parameters can be obtained from file `used_params.txt`

2.1 General

Plugins applied to every sequence, separated by commas. Order is important

1. PluginIndeterminations
2. PluginFindPolyAt
3. PluginAbAdapters
4. PluginUserContaminants
5. PluginContaminants
6. PluginVectors
7. PluginLowQuality
8. PluginLowComplexity
9. PluginExtractInserts

Remove duplicated (clonal) sequences (using CD-HIT 454)

```
remove_clonality: false
```

Minimum insert size for every trimmed sequence

```
min_insert_size_trimmed: 30
```

Minimum insert size for each end of paired-end reads; true paired-ends have both single-ends longer than this value

```
min_insert_size_paired: 40
```

Seqtrim version

```
seqtrim_version: 2.0.67
```

```
min_sequence_size_raw:
```

2.2 Quality

Minimum quality value for every nucleotide

```
min_quality: 20
```

```
window_width:
```

2.3 Contaminants

Blast E-value used as cut-off when searching for contaminations

```
blast_evalue_contaminants: 1.0e-10
```

Minimum required identity (%) for a reliable contamination

```
blast_percent_contaminants: 85
```

Minimum hit size (nt) for considering a true contamination

```
min_contam_seq_presence: 40
```

Genus of input data: contaminations belonging to this genus will be ignored

```
genus:
```

Is a contamination considered a source of sequence rejection? (setting to false will only trim contaminated sequences instead of rejecting the complete read)

```
contaminants_reject: true
Path for contaminants database
    contaminants.fasta
    cont_ribosome.fasta
```

3 Pre-processing statistics

Next figure is equivalent to Figure ?? but using output reads (useful sequences). The mode is expected to decrease but the shape of the plot should be similar.

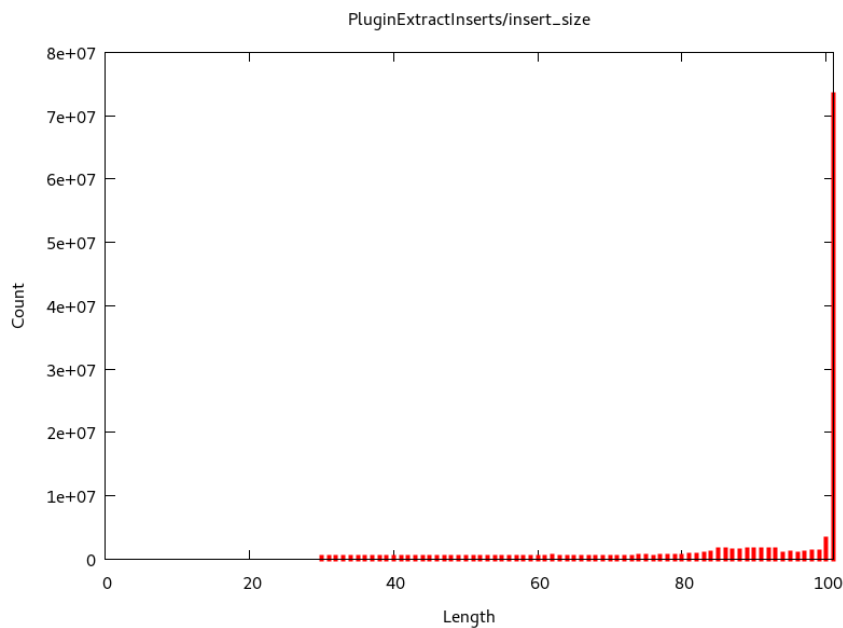


Figure 1: Size distribution of the output sequences. Short sequences ($< \text{min_insert_size_trimmed}$) were removed. [PluginExtractInserts_insert_size.png]

Summary statistics of the SeqTrimNext analysis. Be careful and read all warnings that are indicating concerns about your data. In the files `initial_stats.json` and `stats.json` can be found a full statistics of your data and SeqTrimNext pre-processing

Input reads:	total	133176546
	Smallest read (bp)	101
	Largest read (bp)	101
	Mode (bp)	0
	Mean (bp)	0.0
Output results:	total	5309042
	Rejected	15119179
	Low complexity reads	81487
	Mode (bp)	91
	Mean (bp)	91.9
	Output paired reads	112666838
	Total output reads	117975880
Linkers:		

Table 1: List of the most frequent Vectors found among your reads

Vectors	sequences
Cloning vector pKOHPR7 complete sequence.	569546
Cloning vector pAAV-MCS, complete sequence.	566186
Enterobacteria phage lambda	191613
Cloning vector pVLH/hsp	59190
Cloning vector pWormgate2, complete sequence.	57195

Table 2: List of the most frequent Adapters found among your reads

Adapters	sequences
ABI_Solid3_Adapter_A	379239
ABI_Solid3_GAPDH_Reverse_Primer	327378
ABI_Solid3_GAPDH_Forward_Primer	273922
Illumina_Single_End_Adapter_1	271485
TruSeq_Universal_Adapter	261392

Table 3: List of the most frequent Contaminants found among your reads

Contaminants	sequences
rRNA_small_subunit_Metazoa_Myotis	115601
rRNA_long_subunit_Metazoa_Dasytus	107647
Podospora anserina S mat+ unordered scaffolds, whole genome shotgun sequence	105730
Saccharomyces cerevisiae S288c chromosome XII, complete sequence	102431
Aspergillus niger CBS 513.88 clone An15	74412

Table 4: Summary of nucleotides removed in every plugin.

Plugin	Nucleotides	Percent	Warnings
Low Quality	1825589679	Inf %	OK
Low Complexity	90532333	Inf %	OK
Poly T	39929993	Inf %	OK
Poly A	46766074	Inf %	OK
Contaminants	83515146	Inf %	OK
Adapters	119794561	Inf %	OK
Vectors	70250462	Inf %	OK
Indeterminations	79020	Inf %	OK
Inserts	10854400737	Inf %	iW1

iW1 Warning!, only Inf % of nucleotides are useful

4 Rejected reads

Input sequences	133176546
Output sequences	5309042
Rejected sequences	15119179
Output paired sequences	112666838
Total output sequences	117975880
Low complexity sequences	81487

Table 5: Summary of reads removed in every plugin.

Case	Number of sequences	Percent	Warnings
	15119179	11.353 %	OK
Short inserts	9879950	7.419 %	OK
Empty Inserts	4300675	3.229 %	rdW3
Contaminants	906743	0.681 %	OK
No Valid Inserts	23249	0.017 %	OK
Low Complexity	7647	0.006 %	OK
Indeterminations	810	0.001 %	OK
Unexpected Vector	105	0.000 %	OK
Total rejected	15119179	11.353 %	OK

rdW3 Warning!, a 3.229 % of your sequences are empty (without an insert)

References

- [1] Falgueras et al. SeqTrim: a high-throughput pipeline for preprocessing any type of sequence reads. *BMC Bioinformatics* 11:38 (2010)
- [2] Weizhong Li & Adam Godzik. Cd-hit: a fast program for clustering and comparing large sets of protein or nucleotide sequences. *Bioinformatics* (2006) 22:1658-9

Thanks you for use SeqTrimNext! Send us any comment to scbi support