

SeqTrimNext

Statistics of pre-processing

Plataforma Andaluza de Bioinformática
Universidad de Málaga

April 22, 2020

1 Output Files

SeqTrimNext provides several files, the most interesting ones are in the following directories:

- `output_files`
 - `output.less`, containing an extensive information about the trimming of each sequence. It can be visualised on terminal using the command `less -R`.
 - `used_params.txt`, containing the complete set of parameters used for execution of SeqTrimNext with your data
 - `rejected.txt`, containing a list of rejected sequences together with the reason for their removal.
 - `initial_stats.json`, containing statistics for raw sequences.
 - `stats.json`, containing the statistics of the cleaning process.
 - There is a collection of `folders` that gather sequences with the same MID; each folder contains a `sequences` file (in FASTQ format) with useful reads. There may also exist a file with reads containing low complexity regions. If you want to reconstruct a SFF with the useful segment of each pre-processed read, use `sff_info` file in combination with the original SFF file for the `sfffile` tool.
- `graphs`
 - `size_stats.png`, a graph with the distribution of read lengths in raw data (see Fig. ??).
 - `qualities.png`, a graph to inspect read qualities in raw data (see Fig. ??).
 - `PluginExtractInserts_insert_size.png`, a graph with the distribution of read lengths after SeqTrimNext pre-processing (see Fig. 1).
 - There are other graphs (mostly bar plots) that illustrate the quality of pre-processed reads. All are in PNG format.
- `latex`
 - It is provided as a compressed file `latex.zip` containing all “.tex” files required to compile this document. Graphs are taken from the `graph` folder

2 Relevant parameters

In this section, the relevant parameters used in your experiment are shown. Full information about the parameters can be obtained from file `used_params.txt`

2.1 General

Plugins applied to every sequence, separated by commas. Order is important

1. PluginIndeterminations
2. PluginFindPolyAt
3. PluginAbAdapters
4. PluginUserContaminants
5. PluginContaminants
6. PluginVectors
7. PluginLowQuality
8. PluginLowComplexity
9. PluginExtractInserts

Remove duplicated (clonal) sequences (using CD-HIT 454)

```
remove_clonality: false
```

Minimum insert size for every trimmed sequence

```
min_insert_size_trimmed: 30
```

Minimum insert size for each end of paired-end reads; true paired-ends have both single-ends longer than this value

```
min_insert_size_paired: 40
```

Seqtrim version

```
seqtrim_version: 2.0.67
```

```
min_sequence_size_raw:
```

2.2 Quality

Minimum quality value for every nucleotide

```
min_quality: 20
```

```
window_width:
```

2.3 Contaminants

Blast E-value used as cut-off when searching for contaminations

```
blast_evalue_contaminants: 1.0e-10
```

Minimum required identity (%) for a reliable contamination

```
blast_percent_contaminants: 85
```

Minimum hit size (nt) for considering a true contamination

```
min_contam_seq_presence: 40
```

Genus of input data: contaminations belonging to this genus will be ignored

```
genus:
```

Is a contamination considered a source of sequence rejection? (setting to false will only trim contaminated sequences instead of rejecting the complete read)

```
contaminants_reject: true
Path for contaminants database
    contaminants.fasta
    cont_ribosome.fasta
```

3 Pre-processing statistics

Next figure is equivalent to Figure ?? but using output reads (useful sequences). The mode is expected to decrease but the shape of the plot should be similar.

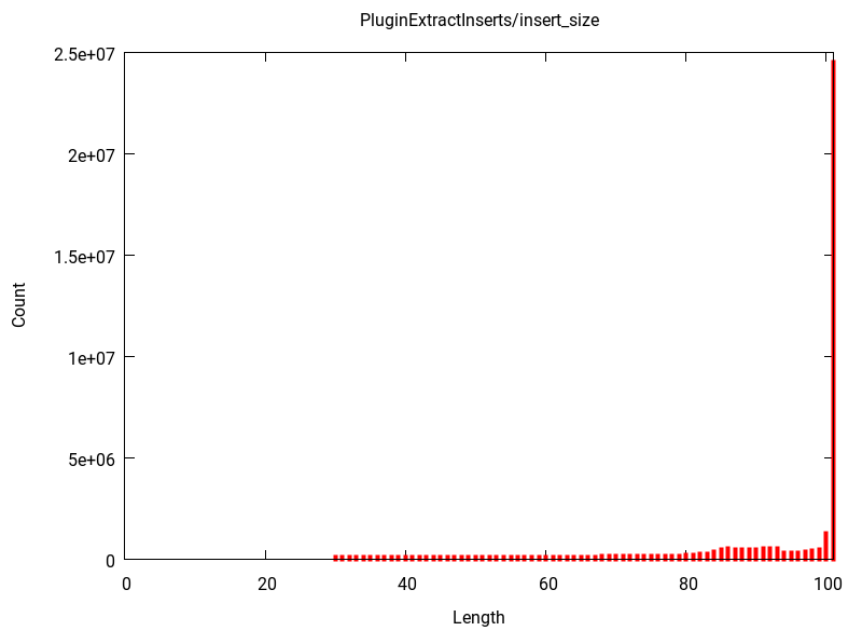


Figure 1: Size distribution of the output sequences. Short sequences ($< \text{min_insert_size_trimmed}$) were removed. [PluginExtractInserts_insert_size.png]

Summary statistics of the SeqTrimNext analysis. Be careful and read all warnings that are indicating concerns about your data. In the files `initial_stats.json` and `stats.json` can be found a full statistics of your data and SeqTrimNext pre-processing

Input reads:	total	45082864
	Smallest read (bp)	101
	Largest read (bp)	101
	Mode (bp)	0
	Mean (bp)	0.0
Output results:	total	1643437
	Rejected	5697087
	Low complexity reads	32680
	Mode (bp)	91
	Mean (bp)	92.5
	Output paired reads	37709660
	Total output reads	39353097
Linkers:		

Table 1: List of the most frequent Vectors found among your reads

Vectors	sequences
Illumina PCR Primer	395560
Illumina Small RNA 3' Adapter	376922
Illumina Multiplexing PCR Primer Index 9	367117
Cloning vector pAAV-MCS, complete sequence.	256524
Cloning vector pKOHPRT complete sequence.	241220

Table 2: List of the most frequent Adapters found among your reads

Adapters	sequences
TruSeq_Adapter_Index_9	337826
TruSeq_Universal_Adapter	290648
TruSeq_Universal_Adapter Illumina_Small_RNA_3p_Adapter_1	252795
Illumina_Small_RNA_3p_Adapter_1	195203
Illumina_Multiplexing_Read2_Sequencing_Primer	131326

Table 3: List of the most frequent Contaminants found among your reads

Contaminants	sequences
rRNA_long_subunit_Metazoa_Dasytus	65837
Saccharomyces cerevisiae S288c chromosome XII, complete sequence	60001
rRNA_small_subunit_Metazoa_Myotis	49808
Podospira anserina S mat+ unordered scaffolds, whole genome shotgun sequence	44154
Schizosaccharomyces pombe 972h- chromosome III, complete sequence	39152

Table 4: Summary of nucleotides removed in every plugin.

Plugin	Nucleotides	Percent	Warnings
Low Quality	561617370	Inf %	OK
Low Complexity	33677611	Inf %	OK
Poly T	12173734	Inf %	OK
Poly A	19605341	Inf %	OK
Contaminants	37937002	Inf %	OK
Adapters	135968001	Inf %	OK
Vectors	60186010	Inf %	OK
Indeterminations	51079	Inf %	OK
Inserts	3644640535	Inf %	iW1

iW1 Warning!, only Inf % of nucleotides are useful

4 Rejected reads

Input sequences	45082864
Output sequences	1643437
Rejected sequences	5697087
Output paired sequences	37709660
Total output sequences	39353097
Low complexity sequences	32680

Table 5: Summary of reads removed in every plugin.

Case	Number of sequences	Percent	Warnings
	5697087	12.637 %	OK
Short inserts	3204731	7.109 %	OK
Empty Inserts	2070007	4.592 %	rdW3
Contaminants	406728	0.902 %	rdW4
No Valid Inserts	6147	0.014 %	OK
Indeterminations	6069	0.013 %	OK
Low Complexity	3340	0.007 %	OK
Unexpected Vector	65	0.000 %	OK
Total rejected	5697087	12.637 %	OK

rdW3 Warning!, a 4.592 % of your sequences are empty (without an insert)

rdW4 Warning!, a 0.902 % of your sequences are from a contaminant organism or from organelles

References

- [1] Falgueras et al. SeqTrim: a high-throughput pipeline for preprocessing any type of sequence reads. *BMC Bioinformatics* 11:38 (2010)
- [2] Weizhong Li & Adam Godzik. Cd-hit: a fast program for clustering and comparing large sets of protein or nucleotide sequences. *Bioinformatics* (2006) 22:1658-9

Thanks you for use SeqTrimNext! Send us any comment to scbi support