

SeqTrimNext

Statistics of pre-processing

Plataforma Andaluza de Bioinformática
Universidad de Málaga

April 22, 2020

1 Output Files

SeqTrimNext provides several files, the most interesting ones are in the following directories:

- `output_files`
 - `output.less`, containing an extensive information about the trimming of each sequence. It can be visualised on terminal using the command `less -R`.
 - `used_params.txt`, containing the complete set of parameters used for execution of SeqTrimNext with your data
 - `rejected.txt`, containing a list of rejected sequences together with the reason for their removal.
 - `initial_stats.json`, containing statistics for raw sequences.
 - `stats.json`, containing the statistics of the cleaning process.
 - There is a collection of `folders` that gather sequences with the same MID; each folder contains a `sequences` file (in FASTQ format) with useful reads. There may also exist a file with reads containing low complexity regions. If you want to reconstruct a SFF with the useful segment of each pre-processed read, use `sff_info` file in combination with the original SFF file for the `sfffile` tool.
- `graphs`
 - `size_stats.png`, a graph with the distribution of read lengths in raw data (see Fig. ??).
 - `qualities.png`, a graph to inspect read qualities in raw data (see Fig. ??).
 - `PluginExtractInserts_insert_size.png`, a graph with the distribution of read lengths after SeqTrimNext pre-processing (see Fig. 1).
 - There are other graphs (mostly bar plots) that illustrate the quality of pre-processed reads. All are in PNG format.
- `latex`
 - It is provided as a compressed file `latex.zip` containing all “.tex” files required to compile this document. Graphs are taken from the `graph` folder

2 Relevant parameters

In this section, the relevant parameters used in your experiment are shown. Full information about the parameters can be obtained from file `used_params.txt`

2.1 General

Plugins applied to every sequence, separated by commas. Order is important

1. PluginIndeterminations
2. PluginFindPolyAt
3. PluginAbAdapters
4. PluginUserContaminants
5. PluginContaminants
6. PluginVectors
7. PluginLowQuality
8. PluginLowComplexity
9. PluginExtractInserts

Remove duplicated (clonal) sequences (using CD-HIT 454)

```
remove_clonality: false
```

Minimum insert size for every trimmed sequence

```
min_insert_size_trimmed: 30
```

Minimum insert size for each end of paired-end reads; true paired-ends have both single-ends longer than this value

```
min_insert_size_paired: 40
```

Seqtrim version

```
seqtrim_version: 2.0.67
```

```
min_sequence_size_raw:
```

2.2 Quality

Minimum quality value for every nucleotide

```
min_quality: 20
```

```
window_width:
```

2.3 Contaminants

Blast E-value used as cut-off when searching for contaminations

```
blast_evalue_contaminants: 1.0e-10
```

Minimum required identity (%) for a reliable contamination

```
blast_percent_contaminants: 85
```

Minimum hit size (nt) for considering a true contamination

```
min_contam_seq_presence: 40
```

Genus of input data: contaminations belonging to this genus will be ignored

```
genus:
```

Is a contamination considered a source of sequence rejection? (setting to false will only trim contaminated sequences instead of rejecting the complete read)

```
contaminants_reject: true
Path for contaminants database
    contaminants.fasta
    cont_ribosome.fasta
```

3 Pre-processing statistics

Next figure is equivalent to Figure ?? but using output reads (useful sequences). The mode is expected to decrease but the shape of the plot should be similar.

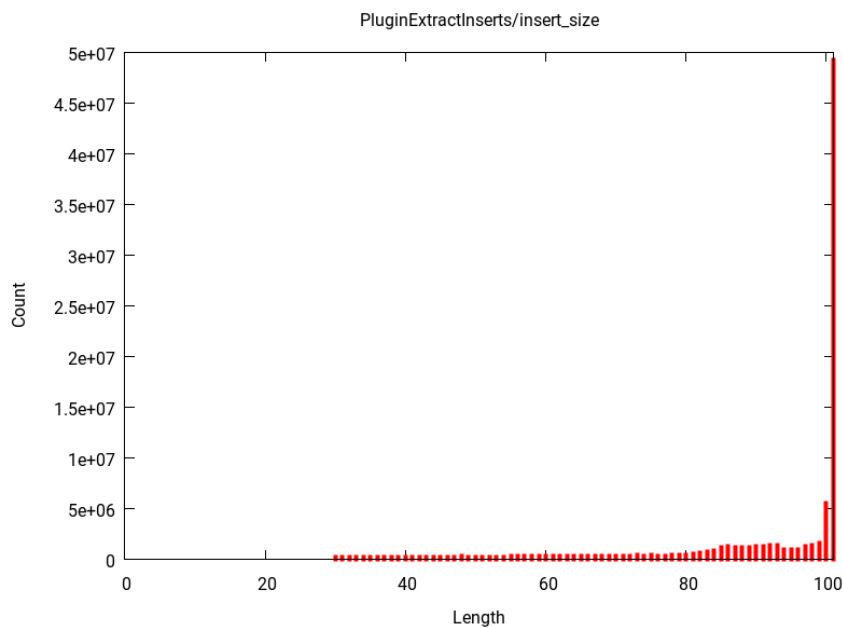


Figure 1: Size distribution of the output sequences. Short sequences ($< \text{min_insert_size_trimmed}$) were removed. [PluginExtractInserts_insert_size.png]

Summary statistics of the SeqTrimNext analysis. Be careful and read all warnings that are indicating concerns about your data. In the files `initial_stats.json` and `stats.json` can be found a full statistics of your data and SeqTrimNext pre-processing

Input reads:	total	101248162
	Smallest read (bp)	101
	Largest read (bp)	101
	Mode (bp)	0
	Mean (bp)	0.0
Output results:	total	3636970
	Rejected	10546342
	Low complexity reads	63764
	Mode (bp)	91
	Mean (bp)	91.5
	Output paired reads	87001086
	Total output reads	90638056
Linkers:		

Table 1: List of the most frequent Vectors found among your reads

Vectors	sequences
Cloning vector pAAV-MCS, complete sequence.	533426
Cloning vector pKOHPRT complete sequence.	522969
Enterobacteria phage lambda	166093
Cloning vector pVLH/hsp	55092
Cloning vector pWormgate2, complete sequence.	46851

Table 2: List of the most frequent Adapters found among your reads

Adapters	sequences
ABISolid3_Adapter_A	288208
TruSeq_Universal_Adapter	286775
ABISolid3_GAPDH_Reverse_Primer	231831
Illumina_Single_End_Adapter_1	200766
ABISolid3_GAPDH_Forward_Primer	197743

Table 3: List of the most frequent Contaminants found among your reads

Contaminants	sequences
rRNA_small_subunit_Metazoa_Myotis	130565
Saccharomyces cerevisiae S288c chromosome XII, complete sequence	105469
rRNA_long_subunit_Metazoa_Dasytus	101010
Schizosaccharomyces pombe 972h- chromosome III, complete sequence	76468
rRNA_small_subunit_Metazoa_Equus	69352

Table 4: Summary of nucleotides removed in every plugin.

Plugin	Nucleotides	Percent	Warnings
Low Quality	1344253540	Inf %	OK
Low Complexity	75988092	Inf %	OK
Poly T	35499173	Inf %	OK
Poly A	37601129	Inf %	OK
Contaminants	70884446	Inf %	OK
Adapters	94385831	Inf %	OK
Vectors	61028255	Inf %	OK
Indeterminations	243990	Inf %	OK
Inserts	8303132392	Inf %	iW1

iW1 Warning!, only Inf % of nucleotides are useful

4 Rejected reads

Input sequences	101248162
Output sequences	3636970
Rejected sequences	10546342
Output paired sequences	87001086
Total output sequences	90638056
Low complexity sequences	63764

Table 5: Summary of reads removed in every plugin.

Case	Number of sequences	Percent	Warnings
	10546342	10.416 %	OK
Short inserts	6907466	6.822 %	OK
Empty Inserts	2853743	2.819 %	rdW3
Contaminants	752038	0.743 %	OK
No Valid Inserts	23668	0.023 %	OK
Low Complexity	6722	0.007 %	OK
Indeterminations	2579	0.003 %	OK
Unexpected Vector	126	0.000 %	OK
Total rejected	10546342	10.416 %	OK

rdW3 Warning!, a 2.819 % of your sequences are empty (without an insert)

References

- [1] Falgueras et al. SeqTrim: a high-throughput pipeline for preprocessing any type of sequence reads. *BMC Bioinformatics* 11:38 (2010)
- [2] Weizhong Li & Adam Godzik. Cd-hit: a fast program for clustering and comparing large sets of protein or nucleotide sequences. *Bioinformatics* (2006) 22:1658-9

Thanks you for use SeqTrimNext! Send us any comment to scbi support